

Phonak ComPilot Air II

User Guide



PHONAK
life is on

Contents

1. Welcome	5
2. Getting to know your ComPilot Air II	6
2.1 Legend	7
3. Getting started	8
3.1 Setting up the power supply	8
3.2 Charging the battery	9
3.3 Switching On or Off	10
3.4 How to wear	10
3.5 Sound demonstration	12
4. Phone calls using ComPilot Air II	14
4.1 How it works	14
4.2 Initial pairing to a phone	15
4.3 Connecting to the phone	17
4.4 Making a call	18
4.5 Accepting a call	19
4.6 Ending a call	20
4.7 Advanced phone options	20
4.8 Additional features	23


5. Streaming audio	25
5.1 Streaming audio via Bluetooth®	25
5.2 Streaming audio from TV	28
6. Using ComPilot Air II as a remote control	29
6.1 Changing the volume	29
6.2 Changing the program	30
6.3 More remote control information	31
7. More detailed information	32
7.1 Main button functionality	32
7.2 Locking the buttons	33
7.3 Switching off the indicator light	34
7.4 Managing multiple inputs	34
7.5 Understanding the indicator lights	36
7.6 Deleting paired Bluetooth devices	38
7.7 Caller identification	39
7.8 VoiceAlerts (spoken messages)	40
7.9 Phonak RemoteControl App	41

8. Troubleshooting	42
9. Service and warranty	50
9.1 Local warranty	50
9.2 International warranty	50
9.3 Warranty limitation	51
10. Compliance information	52
11. Information and explanation of symbols	55
12. Important safety information	59

1. Welcome

Your ComPilot Air II is a premium-quality product developed by Phonak, one of the world's leading companies in hearing technology.

Please read the user instructions thoroughly in order to benefit from all the possibilities your hearing system offers.

 Compatibility information:
Check your hearing aid compatibility with your hearing care professional.

Phonak – life is on

www.phonak-us.com



CE mark applied: 2014

2. Getting to know your ComPilot Air II

The intended use of ComPilot Air II is to receive audio inputs from different Bluetooth® devices and to transmit them wirelessly to your hearing aids.

Virtually any Bluetooth-enabled audio device such as cell phones or MP3 players are supported. Enjoy speech and music from other multimedia devices, such as computers, TVs or music players.




ComPilot Air II can also be used as a remote control to adjust your hearing aids to different listening situations.

i Please note that your hearing care professional must configure your ComPilot Air II before it can be used as a remote control.

Bluetooth® is a registered trademark owned by the Bluetooth SIG, Inc.

2.1 Legend

Buttons

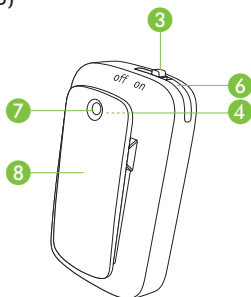
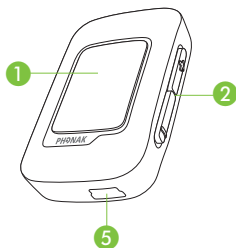
- 1 Main 
- 2 Volume up / down 
- 3 On / Off switch 

Inputs

- 4 Microphone opening
- 5 Charging input (mini-USB)

Others

- 6 Indicator light
- 7 Securing cord eyelet
- 8 Retaining clip



3. Getting started

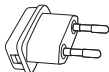
ComPilot Air II must be set up and charged **before** use.

3.1 Setting up the power supply

Select the adapter that corresponds to your country.





US




EU

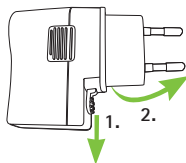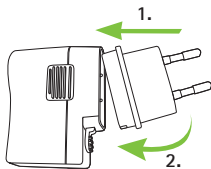- First, insert the top round end of the adapter into the round corner of the power supply.
- Click the bottom end of the adapter in to lock it into place.
- To remove, first pull back the tab on the power supply.
- Pull the adapter upwards to remove it.



3.2 Charging the battery

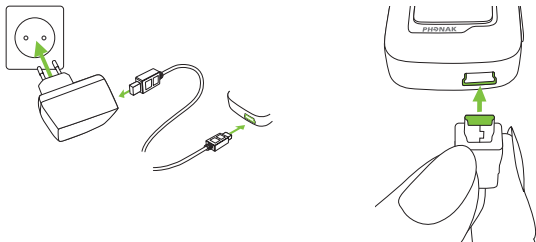- Plug the larger end of the charging cable into the universal power supply. Plug the smaller end into the mini-USB port on ComPilot Air II.
- Plug the power supply into a power outlet.



The indicator light is red during charging and will turn green when charging is complete.


i ComPilot Air II has a non-removable rechargeable battery.

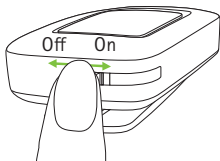
i Before using ComPilot Air II for the first time, charge it for at least 3 hours.

❗ It is safe to leave ComPilot Air II connected to the charger overnight.

❗ The full performance of a new battery is achieved after three complete charge and discharge cycles.

3.3 Switching On and Off

Slide the power switch  as shown to turn ComPilot Air II on or off.



❗ Switch ComPilot Air II off when not in use.

3.4 How to wear

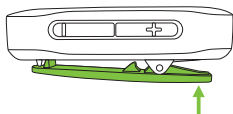
Before any phone call, music or TV streaming:

- Insert your hearing aids into your ears correctly and switch them on.

-
- Switch ComPilot Air II on.



- Clip it securely on to your clothing.



Distance to the mouth:


- The best results are achieved at 20 cm (8").
- A distance of 40 cm (16") to your hearing aids should not be exceeded.

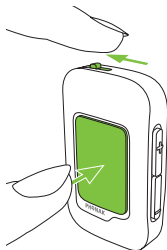


To prevent scratching sounds, do not cover ComPilot Air II with clothing during phone calls.

3.5 Sound demonstration




- ① Wear your hearing aids and clip ComPilot Air II securely to your clothing within 40 cm (16") of your hearing aids (see chapter 3.4).

Press and **hold** the main  button while switching on the device. The indicator light will turn pink.



The hearing aids will switch to streaming mode and you will hear the sound demonstration.



- Use volume   buttons to adjust the volume.
- The demo will stop after 3 minutes or when the main  button is pressed.

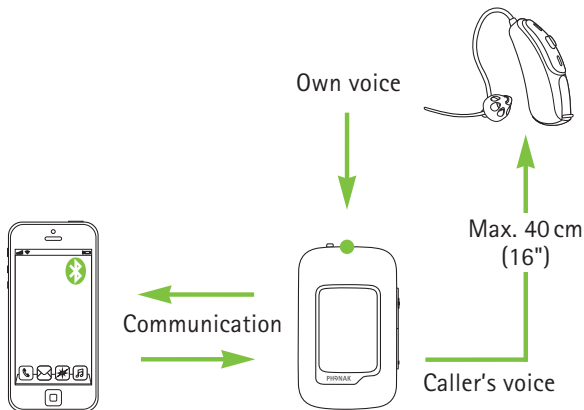
-
- ① Adjusting the volume during the sound demonstration will also set the volume level of the voice alerts (spoken messages). (see chapter 7.8).
 - ① If you cannot hear the sound demonstration, please consult the troubleshooting section.
 - ① The default language setting of ComPilot Air II is English and can be configured by your hearing care professional or the Phonak RemoteControl App (see chapter 7.9).

4. Phone calls using ComPilot Air II

4.1 How it works




ComPilot Air II acts as an interface between your phone and your hearing aids. During a call you will hear the caller's voice directly in your hearing aids.



ComPilot Air II picks up your voice through its own microphone. Speak into ComPilot Air II – not into the phone – when a call is in progress.



4.2 Initial pairing to a phone

i It is only necessary to perform the pairing procedure once with each Bluetooth device. ComPilot Air II will connect automatically to the phone after the initial pairing.

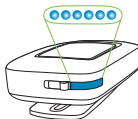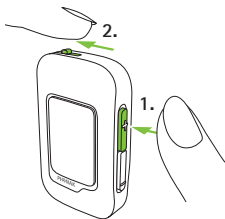- On your phone, ensure that Bluetooth is enabled and search for Bluetooth devices in the connectivity setting menu.
- Switch ComPilot Air II off.
- Press and hold the volume up  button.
- Turn ComPilot Air II on while holding the volume up  button.
- Release the volume up  button when the indicator light starts rapidly blinking blue.



Bluetooth pairing:  + on.

Quick pairing instructions are also shown on the back of ComPilot Air II.

- Start the pairing process on the phone within 1 minute according to the phone user guide.
- Your phone shows a list of found Bluetooth devices. Select "Phonak ComPilot Air II" from this list and enter code "0000" if prompted.
- ComPilot Air II requires access to your phone book. When prompted accept the request.
- After successful pairing, the indicator light will stop rapidly blinking blue.

① For more information about Bluetooth pairing instructions specific to some of the most popular phone manufacturers, go to:
<http://www.phonak-us.com/compilotair>

4.3 Connecting to your phone

Wear your hearing aids and clip ComPilot Air II securely to your clothing within 40 cm (16") distance of your hearing aids and switch them on (see chapter 3.4). After your phone has been paired to ComPilot Air II, they will automatically connect to each other. If they do not connect automatically, go to your phone's Bluetooth menu and select ComPilot Air II from the list.

Please consult your phone's user guide when a manual connection is necessary.

Confirmation of a successful connection is displayed on ComPilot Air II by the indicator light blinking blue.

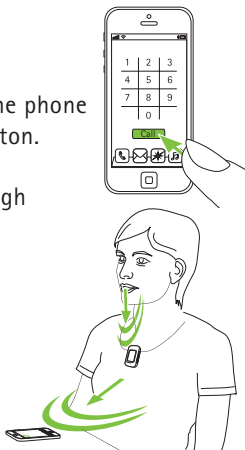
- ❶ The connection will be maintained as long as the devices remain on and within range.
- ❷ ComPilot Air II can connect to two phones at the same time but only one call at a time is possible.
- ❸ When pairing to a phone we recommend disabling the phone key tones and SMS alerts.

4.4 Making a call

Enter the phone number using the phone keypad and press the dialing button.


You will hear the dial tone through your hearing aid.


Speak into ComPilot Air II – not into the phone – when the call is in progress.



4.5 Accepting a call


When a call is received on your phone:

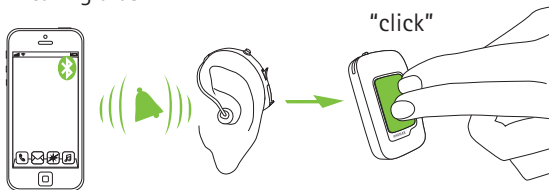- ComPilot Air II indicator light will show a solid blue light.
- It will take about 2 seconds before you hear the ringtone through your hearing aids.
- Press the main  button on your ComPilot Air II once you hear the ringtone through your hearing aids.



If a second call comes in on the same phone:

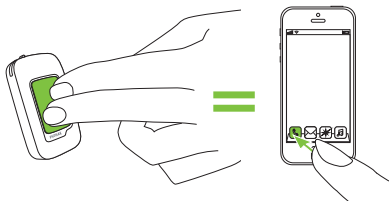
- Press the main  button to accept.
- Press and hold the main  button (2 sec) to reject.

 An incoming call from a second paired phone will be ignored during an active call.

4.6 Ending a call

Press the main  button to end the call.


You can also end a call on your phone instead of your ComPilot Air II.



You will hear the previously active program or audio source.

4.7 Advanced phone options

Rejecting a call



Press and hold the main  button (2 sec) until the ringtone stops. Using your phone to reject the call will have the same effect.

Setting your phone volume

Use your phone's volume buttons to adjust the phone volume to a comfortable level.

Surrounding sounds picked up by your hearing aids are reduced according to the settings made by your hearing care professional.

For an additional reduction of surrounding sounds:

- Press and **hold** the volume down  button 2 sec for more reduction.
- Press and **hold** the volume up  button 2 sec to return to the default reduction.

Disabling unwanted alerts

Some phones are configured to play key-press tones, SMS announcements and confirmation beeps.

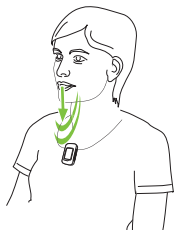
You are able to hear these alerts in your hearing aids.

If this is distracting, refer to your phone user guide for information about disabling these alerts.

How to wear ComPilot Air II during phone calls

- During phone calls, wear ComPilot Air II clipped securely to your clothing as shown.

Note: Maximum distance between ComPilot Air II and the phone is 5–10 m (15–30 ft).



- In very noisy environments, you can move ComPilot Air II closer to your mouth.



4.8 Additional features

The following features can be set either using the Phonak RemoteControl App (see chapter 7.9) or by your hearing care professional.

i Your phone must support these features. Please consult your phone user manual for details.

Caller's name announcement


By default, ComPilot Air II will announce the caller's name as shown on the phone display (caller identification).

The following features can be enabled separately.

Press and hold the Main button (2 sec.) to use them:

1. *Last number redial*: To redial the last dialed number.
2. *Voice dial*: Allows you to say the person's name instead of typing the number.
3. *Dial preferred number*: Directly dials a predefined number.

-
4. *Hold call:* While in a phone conversation, this will put the call on hold. The caller stays on the line but there is no audio transfer and the hearing aids switch to the previously active hearing program.
 5. *Transfer call:* This allows you to hand your phone to another person. Your hearing aids will switch back to the previously active hearing program.

To transfer the call back to ComPilot Air II, press and hold the main  button. The call will now resume through your hearing aids.

5. Streaming audio

- i** Before streaming, wear your hearing aids and ComPilot Air II correctly and switch them on (see 3.4).

5.1 Streaming audio via Bluetooth



Streaming audio via cell phone

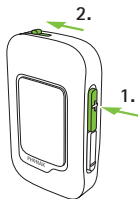
You can stream music wirelessly from your cell phone:


- Make sure ComPilot Air II is connected to your phone as described in chapter 4.3.
- Use your phone music player to play music.

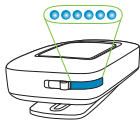
Streaming audio via other devices

To use ComPilot Air II with a Bluetooth enabled audio device (e.g. MP3 player), you must first pair the devices:

- Switch the ComPilot Air II off.
- Press and hold the volume up  button.
- Switch ComPilot Air II on while holding the volume up  button.




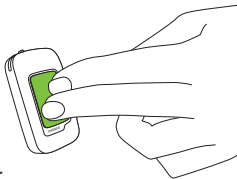
-
- Release the volume up  button when the indicator light starts rapidly blinking blue.
 - Start the pairing process on the Bluetooth device within 1 minute according to the device's user guide and use code "0000" if prompted.
 - After successful pairing, the indicator light will stop rapidly blinking.
 - Press play on the Bluetooth device if needed.
 - Streaming will start automatically and the indicator light will turn solid blue displaying the active streaming.




To adjust the volume of the audio signal use the volume controls on the audio device.

Press the main  button to pause and restart streaming.

To stop Bluetooth streaming completely, press and **hold** the main  button (2 sec), switch the other Bluetooth device off or disable it's Bluetooth feature.



-
- ① Please refer to the device's user guide or dealer support if you cannot complete the pairing procedure with the generic steps described.
 - ① Maximum ComPilot Air II streaming time is 4 hours with a fully charged battery.

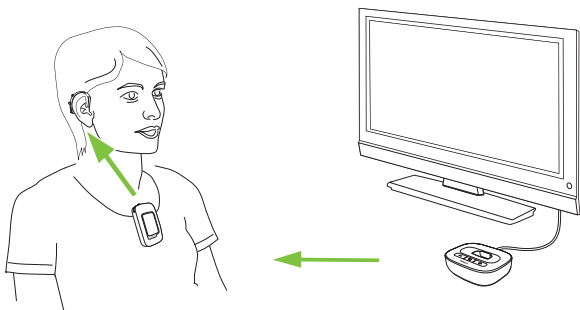
Note: During streaming, ComPilot Air II cannot be used as a hearing aid remote control. Press and hold the main  button (2 sec) to stop the Bluetooth stream completely.

5.2 Streaming audio from TV

It is recommended to use the dedicated and optimized Phonak TVLink II basestation for watching TV or listening to music through a television.

It offers the following benefits:

- Wireless transmission up to 30 meters (90 ft).
- Analog and digital inputs (optical, coaxial).
- Charging cradle for ComPilot Air II.



i Refer to the TVLink II user guide for pairing instructions.

6. Using ComPilot Air II as a remote control

- ❗ Before this functionality can be used, your hearing care professional must enable the remote control in ComPilot Air II.



If ComPilot Air II is used as a handheld remote control:

- Max. distance to the hearing aids is 55 cm / 22".
- Hold ComPilot Air II as shown.



6.1 Changing the volume

Press buttons:



-  to increase the volume.
-  to decrease the volume.




- ❗ Adjusting the volume with ComPilot Air II will have the same effect as adjusting the volume directly on your hearing aids.


Reducing surrounding sounds

You can reduce the surrounding sounds with a single button press. This may be helpful during music or phone streaming. The level of reduction is set by your hearing care professional.

- Press and hold the volume down  button (>2 sec) for more reduction.
- Press and hold the volume up  button (>2 sec) to return to the original hearing aid volume.


6.2 Changing the program

 Changing hearing programs is only possible if no streaming is in progress.

- Press the main  button to proceed to the next available program.




6.3 More remote control information

- Press and hold the main  button (>2 sec) to select your favorite hearing program or jump to your start-up program configured by your hearing care professional or using the Phonak RemoteControl App.
- Both hearing aids are controlled at the same time.

7. More detailed information

7.1 Main button functionality

The Main  button has multiple functionalities:



1. *Phone*: It accepts and ends a phone call.
2. *Stream*: It pauses and resumes the streaming from an audio device. A long press (2 sec) will disconnect the streaming and programs can be changed again.
3. *Remote control*: Used to change the hearing programs.

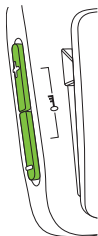


- ❗ The hearing programs cannot be changed when:
- a phone call is in progress
 - a Bluetooth device is streaming or paused
 - the remote control feature has not been enabled by your hearing care professional

7.2 Locking the buttons

To prevent the buttons from being accidentally pressed:

- Press the volume up  and down  buttons simultaneously to lock or unlock the buttons.
- The indicator light will blink red 3 times, indicating that the buttons are locked.
- When locked, the indicator light will blink red 3 times and you may hear three beeps through your hearing aids each time a button is pressed.



i Accepting an incoming call is still possible.

i In locked state, the indicator light will stop blinking.

i The key symbol on the back reminds you of the lock function.

7.3 Switching off the indicator light

If the indicator light is distracting:

- Lock the ComPilot Air II as described in chapter 7.2.
- During the locked state the indicator light remains switched off.

i A low battery indicator will still be shown.

7.4 Managing multiple inputs




When an audio source is connected via ComPilot Air II, your hearing aids will automatically switch to the special streaming programs.

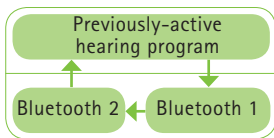
It is possible to have two Bluetooth audio sources connected.

i A phone call always has priority.
You will hear the ringtone even when listening to other audio sources. If you reject or terminate the call, audio streaming will resume.

i If you want to listen to music from your connected cell phone, you must start the music via your phone initially; it cannot be started from ComPilot Air II.







Selecting audio sources

- *One source:* Each press of the main  button will pause or resume the source.
- Press and hold the main  button (2 sec) to completely stop the music streaming. It must be re-started via the phone menu.
- *Multiple sources:* After each press of the main  button, ComPilot Air II selects the next connected source in a fixed sequence.
- Every connected audio source has a fixed place in the sequence, as shown in the picture.



7.5 Understanding the indicator light



Bluetooth pairing	Rapid blink blue 
1 Bluetooth device connected	Short blink blue 
2 Bluetooth devices connected	Double-blink blue 
Phone or music via Bluetooth	Solid blue 
Demo sound playing	Solid pink 
Charging	Solid Red 

Charged and full battery

Solid green



Switching On

Green 2 seconds



Switching Off

Solid red for 2 seconds



Battery above 10%

Slow blink green



Less than 10% battery left.
Please charge

Short flash red



Enable /
disable lock state

Triple flash red (locked)






Triple flash green (unlocked)



7.6 Deleting paired Bluetooth devices

If you encounter issues with one of the paired devices, you may consider deleting the pairings from ComPilot Air II.

 Deleting the pairing of Bluetooth devices removes the ability to connect to any previously paired Bluetooth device.

1. Turn ComPilot Air II on while pressing the volume up  button.
2. Release the volume up  button only after 10 seconds when the rapidly blue blinking indicator light turns green.
3. Delete the ComPilot Air II pairing entries on your mobile phone before any re-pairing.
4. Repeat the pairing procedure as described in chapter 4.2.

7.7 Caller identification

ComPilot Air II is able to tell you the name of the caller.

ComPilot Air II will try to access information about the caller from your phone book. It will be spoken through your hearing aids as follows:

- The stored name in your phone's phonebook.
- The number if the name is not known.

Your hearing care professional can configure the language setting of your ComPilot Air II.

You can also use the Phonak RemoteControl App to configure the language setting (see chapter 7.9).

7.8 VoiceAlerts (spoken messages)

VoiceAlerts inform you about the different operation modes and status of ComPilot Air II using spoken messages.

Your hearing care professional can configure the language setting as well as the speech levels:

- *Warnings only:* Spoken information about a low ComPilot Air II battery and during the Bluetooth pairing process.
- *Standard:* Additional information about switching ComPilot Air II on or off, start or end charging during streaming, input source selection announcements.
- *Detailed:* Additional information about start or end of charging even if no streaming source is active, Bluetooth device connect or disconnect and transfer or hold a call.

7.9 Phonak RemoteControl App

The Phonak RemoteControl App is a Smartphone application, which offers additional control of your hearing aids.



Features are shown on your Smartphone screen and give you the ability to configure your hearing system.

Experience more:

- Individual left / right volume control.
- Direct program or alternative audio input selection.
- Support for the Bluetooth pairing process.

① For more information and to download the RemoteControl App visit the website:
<http://www.phonak-us.com/rcapp>



8. Troubleshooting

For any problems not listed in this user guide, please contact your hearing care professional.

Causes

What to do

During Bluetooth pairing my other device is prompting for a code, what do I have to enter?

- | | |
|---|---|
| <ul style="list-style-type: none">■ The pairing code is required for the initial pairing of ComPilot Air II with your Bluetooth device. | <ul style="list-style-type: none">■ Enter "0000" (four zeros) when prompted for a code. |
|---|---|

I hear intermittency during an audio transfer via ComPilot Air II.

- | | |
|---|--|
| <ul style="list-style-type: none">■ The Bluetooth device is too far away, or worn in a back pocket (body shielding).■ ComPilot Air II searches for known devices for 2 minutes after switching on.■ A phone connected to other devices such as a PC may not be able to continuously stream music. | <ul style="list-style-type: none">■ Do not exceed the maximum distance 5–10 m (15–30 ft.) and keep Bluetooth devices in front of you.■ After 2 minutes, ComPilot Air II will stop searching and the intermittency will stop.■ When using your Bluetooth phone as music player, disconnect it from other devices such as PCs. |
|---|--|

Causes

What to do


The volume of the music or phone conversation is too loud or too quiet.

- Different audio sources have different outputs.

- To change the volume of the phone, or any other sound source to ComPilot Air II, use the volume control of that external device.


I have difficulties understanding the phone conversation in a noisy environment.

- Your hearing aid microphones are set too high during the call.

- Reduce the environmental noise by pressing the  button for 2 seconds (see chapter 6.1).
- Use the phone volume button to increase the phone volume.

The caller can hear me but I cannot hear the caller.

- The call was transferred to your phone.
- Some phones, when accepting the call through the phone button, don't use ComPilot Air II.

- Make sure the call is transferred back to ComPilot Air II via the phone menu.
- Always accept the call by pressing the main  button.

Causes

What to do

ComPilot Air II is no longer recognized by my phone or other Bluetooth devices to which it was previously paired.

- | | |
|--|---|
| <ul style="list-style-type: none">■ 2 minutes after not finding a device, ComPilot Air II will stop seeking devices to save power.■ The devices are too far away from each other.■ ComPilot Air II can be paired to a maximum of 6 different Bluetooth devices. If the internal memory is full, new pairings may overwrite previous pairings.■ Pairing was deleted. | <ul style="list-style-type: none">■ Switch ComPilot Air II off. Wait until the red indicator light switches off, then turn on again.■ Bring devices within 1 meter (3 ft) range and try again.■ Repeat the pairing between ComPilot Air II and the Bluetooth device which was overwritten.■ Repeat the pairing between ComPilot Air II and the Bluetooth device which was deleted. |
|--|---|

I am not streaming but my hearing aids keep switching between the streaming and the normal hearing aid program.

- | | |
|---|--|
| <ul style="list-style-type: none">■ Your phone sends sounds to ComPilot Air II due to incoming SMS or mails.■ Phone key tones are switched on. | <ul style="list-style-type: none">■ Disable the system sounds for reminders, alarms, SMS, etc., via your phone menu.■ Switch off phone key tones via your phone menu. |
|---|--|


Causes

What to do

I hear my phone's key-press tones in my hearing aids.

- | | |
|---|---|
| ■ The key-press tones may be enabled on your phone. | ■ Deactivate key-press and confirmation tones on your phone. Please refer to your phone's user guide. |
|---|---|

Every time I switch on my phone, the music player starts.

- | | |
|---|--|
| ■ For some phones this is normal behavior and not initiated by ComPilot Air II. | ■ Some phones can be configured not to start the music player automatically.
■ Press and hold the main  button to stop streaming.
■ Stop the music player via your phone menu after the phone is connected. |
|---|--|

Every time I get in my car the ComPilot Air II connection to my phone behaves unexpectedly.

- | | |
|--|--|
| ■ The car Bluetooth system connects to your phone. | ■ Consider disconnecting the car system from your phone. |
|--|--|

Causes

What to do

My phone is ringing but the indicator light is not solid blue and there is no ringing signal through my hearing aids.

- | | |
|---|---|
| <ul style="list-style-type: none">■ Bluetooth is disabled on your phone.■ ComPilot Air II and your phone are not connected to each other.■ ComPilot Air II is not paired to your phone.■ The distance is too great between ComPilot Air II and your phone. | <ul style="list-style-type: none">■ Ensure Bluetooth is enabled according to your phone's user guide.■ 1. Disable the Bluetooth functionality on your phone and turn it on again.
2. Restart ComPilot Air II by switching it off and then on. The Bluetooth functionality will be automatically reactivated.
3. Activate ComPilot Air II by selecting it in the phone's "Active device" menu.
4. Reduce the distance to phone.■ Follow the pairing process in this guide.■ Wear ComPilot Air II clipped to your clothing. Keep distance to phone <10 m (<30 ft). |
|---|---|

Causes

What to do

My phone is ringing, the indicator light is solid blue but there is no ringing signal through my hearing aids.

- | | |
|--|--|
| <ul style="list-style-type: none">■ ComPilot Air II is out of range of your hearing aids.■ Phone is set to vibrate.■ You are already on a phone call and a second call cannot be recognized. | <ul style="list-style-type: none">■ Wear ComPilot Air II correctly clipped to your clothing with max. 40 cm (16") distance to your hearing aids.■ Enable ring tone on phone.■ Hang up the first phone call, then accept the second call. |
|--|--|

My phone rings but the name of the caller is not heard acoustically.

- | | |
|---|---|
| <ul style="list-style-type: none">■ The caller identification feature is not available or not enabled on your phone.■ Two phones are connected to your ComPilot Air II.■ Caller identification may be switched off. | <ul style="list-style-type: none">■ Please refer to your phone's user guide to check the availability of this feature.■ If two phones are connected caller identification is only available on the phone that was last paired to ComPilot Air II.■ Please ask your hearing care professional to enable caller identification on your ComPilot Air II. |
|---|---|

Causes

What to do

During a phone call the other party has difficulty hearing me.

- | | |
|--|--|
| <ul style="list-style-type: none">■ ComPilot Air II microphone opening may be covered.■ ComPilot Air II may not be worn correctly.■ ComPilot Air II may rub against clothing.■ The surrounding noise may be too loud. | <ul style="list-style-type: none">■ Make sure the microphone opening is not covered by part of your body, clothing or dirt and debris.■ Make sure your ComPilot Air II is clipped securely to your clothing within 40 cm (16") distance of your hearing aids.■ Reduce movement while you speak.■ Bring ComPilot Air II closer to your mouth.■ Suggest to your callers that they increase their phone volume. |
|--|--|

ComPilot Air II is unresponsive or there is no indicator light when switching On.

- | | |
|--|--|
| <ul style="list-style-type: none">■ ComPilot Air II battery may be dead.■ It may indicate a software problem. | <ul style="list-style-type: none">■ Charge for at least 1 hour.■ Switch ComPilot Air II off, wait 30 seconds, then switch On again. |
|--|--|

Causes

What to do

The operating time of my ComPilot Air II is reduced substantially.

- | | |
|--|--|
| <ul style="list-style-type: none">■ The battery typically needs several charging cycles to achieve full performance. | <ul style="list-style-type: none">■ Charge for at least 3 hours the first time you charge ComPilot Air II.■ Let a new ComPilot Air II discharge completely three times. |
|--|--|

I have interruptions while my audio streaming is active.

- | | |
|---|---|
| <ul style="list-style-type: none">■ ComPilot Air II is too far away from your hearing aids while streaming.■ Bluetooth signal is out of range. | <ul style="list-style-type: none">■ Wear ComPilot Air II clipped to your clothing with max. 40cm (16") distance of your hearing aids.■ Keep the Bluetooth device within the operating range. |
|---|---|

I cannot hear the demo sound.

- | | |
|--|---|
| <ul style="list-style-type: none">■ Hearing aids are not worn correctly.■ ComPilot Air II is not worn correctly or is switched off. | <ul style="list-style-type: none">■ Make sure your hearing aids are inserted correctly.■ Wear ComPilot Air II clipped to your clothing with max. 40cm (16") distance of your hearing aids.■ Make sure ComPilot Air II is charged and switched on. |
|--|---|

9. Service and Warranty

9.1 Local warranty

Please ask your hearing care professional about the terms of the local warranty.

9.2 International warranty

Phonak offers a one year limited international warranty, valid as of the date of purchase. This limited warranty covers manufacturing and material defects. The warranty only comes into force if proof of purchase is shown. The international warranty does not affect any legal rights that you might have under applicable national legislation governing sale of consumer goods.

9.3 Warranty limitation

This warranty does not cover damage from improper handling or care, exposure to chemicals, immersion in water or undue stress. Damage caused by third parties or non-authorized service centers renders the warranty null and void. This warranty does not include any services performed by a hearing care professional in their office.

Serial number:

Date of purchase:


Authorized hearing
care professional
(stamp / signature):

10. Compliance information

Declaration of conformity

Hereby Phonak AG declares that this Phonak product is in compliance with the essential requirements of the Medical Devices Directive 93/42/EEC as well as the Radio and Telecommunications Terminal Equipment Directive 1999/5/EC. The full text of the Declaration of Conformity can be obtained from the manufacturer or the local Phonak representative whose address can be taken from the list on <http://www.phonak-us.com> (worldwide locations).

Australia:

Supplier Code Number	 N15398
----------------------	--

New Zealand:

Supplier Code Number	Z1285
----------------------	-------

Notice 1:

This device complies with Part 15 of the FCC Rules and with RSS-210 of Industry Canada. Operation is subject to the following two conditions:

- 1) this device may not cause harmful interference, and
- 2) this device must accept any interference received, including interference that may cause undesired operation.

Notice 2:

Changes or modifications made to this device not expressly approved by Phonak may void the FCC authorization to operate this device.

Notice 3:

This device has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules and ICES-003 of Industry Canada. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This device generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this device does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the device and receiver.
- Connect the device into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio / TV technician for help.

Notice 4:

Japanese Radio Law and Japanese Telecommunications Business Law Compliance. This device is granted pursuant to the Japanese Radio Law (電波法) and the Japanese Telecommunications Business Law (電気通信事業法). This device should not be modified (otherwise the granted designation number will become invalid)."

11. Information and explanation of symbols



With the CE symbol, Phonak AG confirms that this Phonak product meets the requirements of the Medical Devices Directive 93/42/EEC as well as the R&TTE Directive 1999/5/EC on radio and telecommunications equipment.



This symbol indicates that the products described in these user instructions adhere to the requirements for an applied part of Type B of EN 60601-1. The surface of the device is specified as an applied part of Type B.



This symbol indicates that it is important for the user to read and take into account the relevant information in this user guide.



Indicates the medical device manufacturer, as defined in EU Directive 93/42/EEC.



This symbol indicates that it is important for the user to pay attention to the relevant warning notices in these user guides.



Important information for handling and product safety.



EMC and radio communications compliance label Australia.

Operating conditions

The device is designed such that it functions without problems or restrictions if used as intended, unless otherwise noted in this user guide.



Temperature transport and storage:
-20°C to +60°C
(-4° to +140° Fahrenheit)
Temperature operation:
0° to 40° Celsius
(+32° to +104° Fahrenheit)



Keep dry.



Humidity transportation and storage:
<90% (non condensing).
Humidity operation: <95% (non
condensing).



Atmospheric pressure: 200 hPA to
1500 hPA.



The symbol with the crossed-out garbage bin is to make you aware that this device may not be thrown away as normal household waste. Please dispose of old or unused device, at waste disposal sites intended for electronic waste, or give your device to your hearing care professional for disposal. Proper disposal protects the environment and health.



The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Phonak is under license. Other trademarks and trade names are those of their respective owners.





Japanese mark for certified radio equipment.


12. Important safety information

The following pages contain important safety information.

Hazard warnings

 Keep this device out of reach of children under 3 years.

 Phonak devices may generate magnetic fields. Tests that Phonak conducted showed that no interferences with the Implantable devices (e.g. pacemakers, defibrillators, etc.) were caused that would impact the safe and effective use of the Implantable devices. Users of Implantable devices should contact their doctor and/or the manufacturer of the Implantable device before they start using a Phonak device. Stop using the Phonak device if any influence of the Phonak device on the Implantable device is experienced and contact the manufacturer of the Implantable device for advice.

 Only use hearing aids that have been specially programmed for you by your hearing care professional.

⚠ Changes or modifications to the device that were not explicitly approved by Phonak are not permitted.

⚠ Opening the device might damage it. If problems occur which cannot be resolved by following the remedy guidelines in the troubleshooting section of this user guide, consult your hearing care professional.

⚠ This product has an embedded, non-replaceable battery. Do not attempt to open the product or remove the battery as this may cause injury and damage the product.

⚠ Dispose of electrical components in accordance with your local regulations by Phonak AG.

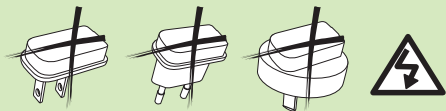
⚠ External devices may only be connected if they have been tested in accordance with corresponding IECXXXXX standards. Only use accessories approved by Phonak AG, to avoid any electrical shock.

⚠ Using the device, its connectors or cables in any way contradictory to their intended purpose as described in this user guide (e.g., wearing the USB cable around the neck) can cause injury.

⚠ When operating machinery, ensure that no parts get caught in the machine.

⚠ The mini-USB port is to be used for the described purpose only.

⚠ For safety reasons, only use chargers supplied by Phonak or stabilized chargers with a rating of 5 VDC, min. 500 mA.



⚠ Caution: electric shock. Do not insert plug alone into electrical outlets.

⚠ Do not use the device in explosive areas (mines or industrial areas with danger of explosions, oxygen-rich environments or areas where flammable anesthetics are handled) or where electronic equipment is prohibited.

Information on product safety

- ① Protect the device from excessive moisture (bathing, swimming), heat (radiator, car dashboard) or direct skin contact when sweating (workout, fitness, sport).
- ① Do not drop the device. Dropping onto a hard surface can damage the device.
- ① Special medical or dental examination including radiation described below, may adversely affect the correct functioning of your device. Remove and keep it outside the examination room/area before undergoing:
 - Medical or dental examination with X-ray (also CT scan).
 - Medical examinations with MRI/NMRI scans, generating magnetic fields.
- ① Protect the device, charging slot and power supply from dirt and debris.

-
- ① Clean the device using a damp cloth. Never use household cleaning products (detergent, soap, etc.) or alcohol to clean the device. Never use a microwave or other heating devices to dry the device.
 - ① The digitally coded, inductive transmission technology used in the device is highly reliable and experiences virtually no interference from other devices. It should be noted, however, that when operating the hearing system near computer equipment, larger electronic installations or other strong electromagnetic fields, it may be necessary to be at least 60 cm (24") away from the interfering device to ensure proper operation.
 - ① Keep the device at least 10 cm away from any kind of magnets.
 - ① Do not connect an USB cable exceeding 3 meters (9 ft) in length to the device.
 - ① When the device is not in use, turn it off and store it safely.

Other important information

- ① High-powered electronic equipment, larger electronic installations and metallic structures may impair and significantly reduce the operating range.
- ① If the hearing aids do not respond to the device because of an unusual field disturbance, move away from the disturbing field.

[illegible]



Phonak LLC
4520 Weaver Parkway
Warrenville, IL 60555

www.phonak-us.com

