

A Phonak is for everyone






Solutions for Better Hearing – New Products

Spring 2010



Versáta Art

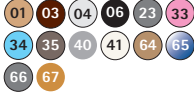

	Features	CIC/MC	ITC/HS	FS
		CIC/MC Power	ITC/HS Power	FS Power
		Petite	Petite/Wireless	Wireless
Fitting Recommendations	Fitting Range	35 dB/113 dB SPL - CIC/MC 55 dB/115 dB SPL - CIC/MC P	48 dB/113 dB SPL - ITC/HS 55 dB/115 dB SPL - ITC/HS P	50 dB/115 dB SPL - FS 55 dB/122 dB SPL - FS P
<ul style="list-style-type: none"> The best solution for moderately active lifestyles Styles for every preference from CIC to micro open fits Connectivity to devices such as cell phones, iPod™, MP3, TV or GPS systems 	Battery Size	10A	312	13
	SoundRecover <i>non-linear frequency compression</i>	■	■	■
	Programmable Channels	16	16	16
	SoundFlow Advanced <i>adapts for seamless transition</i>	■	■	■
	VoiceZoom <i>automatic adaptive up to 12 noise sources</i>		■	■
Standard Instrument Specific Features	RealEar Sound - <i>restores sound localization to determine source of sound</i>			
<ul style="list-style-type: none"> SelfLearning and DataLogging with SmartTips (<i>records and stores personal use and listening preferences</i>) End of battery life indicator 	WhistleBlock <i>more stable amplification without sound quality impact</i>	■	■	■
	BassBoost <i>increased low frequency gain and output</i>	■ <i>(CIC/MC P only)</i>	■ <i>(ITC/HS P only)</i>	■ <i>(FS P only)</i>
	WindBlock <i>suppresses wind noise without microphone changes</i>	■	■	■
	NoiseBlock <i>detects and reduces annoying sounds</i>	■ <i>(16 Channels)</i>	■ <i>(16 Channels)</i>	■ <i>(16 Channels)</i>
Options	AOV - Acoustically Optimized Vent <i>minimizes occlusion with maximum vent size</i>	■	■	■
<ul style="list-style-type: none"> Choice of BTE colors:  Choice of micro BTE colors:  Choice of micro Petite BTE colors:  Choice of Custom colors: See page 8 	EasyPhone - <i>Automatic Telephone</i>		■ <i>(Optional)</i>	■ <i>(Optional)</i>
	Telecoil		■	■
	Number of Manual Programs	4	4	4
	EasyAudio - <i>binaural stereo streaming</i> EasyBluetooth™ - <i>handsfree cell phone use</i>		■ <i>(Wireless only)</i>	■
	QuickSync <i>binaural volume and program control</i>		■ <i>(Wireless only)</i>	■
	iCom - <i>wireless Bluetooth™ & FM</i>		■ <i>(Wireless only)</i>	■
	iCube <i>CableFree Fitting</i>		■ <i>(Wireless only)</i>	■
	myPilot (remote control) <i>Bi-directional access</i>		■ <i>(Wireless only)</i>	■
	KeyPilot2 & WatchPilot2 (remote controls)	■	■	■
	Design Integrated FM Receiver			
	Universal FM Receiver <i>myLink</i>		■	■
	FM Transmitters <i>inspiro, SmartLink+, EasyLink+, ZoomLink+</i>		■	■
	Tamper proof Battery compartment			


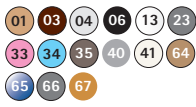
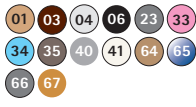
Versáta Art



Features	BTE M	BTE P	BTE SP	Micro
	Wireless	Wireless	Wireless	Wireless/Petite
Fitting Range	53 dB/125 dB SPL	63 dB/130 dB SPL	70 dB/131 dB SPL	46 dB/123 dB SPL - Petite 53 dB/124 dB SPL - Wireless
Battery Size	13	13	13	10A - Petite 13 - Wireless
SoundRecover <i>non-linear frequency compression</i>	■	■	■	■
Programmable Channels	16	16	16	16
SoundFlow Advanced <i>adapts for seamless transition</i>	■	■	■	■
VoiceZoom <i>automatic adaptive up to 12 noise sources</i>	■	■	■	■
RealEar Sound - <i>restores sound localization to determine source of sound</i>	■	■	■	■
WhistleBlock <i>more stable amplification without sound quality impact</i>	■	■	■	■
BassBoost <i>increased low frequency gain and output</i>			■	
WindBlock <i>suppresses wind noise without microphone changes</i>	■	■	■	■
NoiseBlock <i>detects and reduces annoying sounds</i>	■ (16 Channels)	■ (16 Channels)	■ (16 Channels)	■ (16 Channels)
AOV - Acoustically Optimized Vent <i>minimizes occlusion with maximum vent size</i>				
EasyPhone - <i>Automatic Telephone</i> EasyFM - <i>Automatic FM Activation</i>	■	■	■	■ (Wireless only)
Telecoil	■	■	■	■
Number of Manual Programs	4	4	4	4
EasyAudio - <i>binaural stereo streaming</i> EasyBluetooth™ - <i>handsfree cell phone use</i>	■	■	■	■ (Wireless only)
QuickSync <i>binaural volume and program control</i>	■	■	■	■ (Wireless only)
iCom - <i>wireless Bluetooth™ & FM</i>	■	■	■	■ (Wireless only)
iCube <i>CableFree Fitting</i>	■	■	■	■ (Wireless only)
myPilot (remote control) <i>Bi-directional access</i>	■	■	■	■ (Wireless only)
KeyPilot2 & WatchPilot2 (remote controls)	■	■	■	■
Design Integrated FM Receiver	ML9i	ML9i	ML9i	ML12i (Wireless only)
Universal FM Receiver <i>myLink</i>	■	■	■	■
FM Transmitters <i>inspiro, SmartLink+, EasyLink+, ZoomLink+</i>	■	■	■	■
Tamper proof Battery compartment	■ (Optional)	■ (Optional)	■ (Optional)	■ (Wireless option only)

Certéna Art



	Features	CIC/MC CIC/MC Power	ITC/HS ITC/HS Power	FS FS Power
		Petite	Petite/Wireless	Wireless
Fitting Recommendations <ul style="list-style-type: none"> The ideal solution for calm lifestyles Styles for every preference from CIC to micro open fits Connectivity to devices such as cell phones, iPod™, MP3, TV or GPS systems 	Fitting Range	35 dB/113 dB SPL - CIC/MC 55 dB/115 dB SPL - CIC/MC P	48 dB/113 dB SPL - ITC/HS 55 dB/115 dB SPL - ITC/HS P	50 dB/115 dB SPL - FS 55 dB/122 dB SPL - FS P
	Battery Size	10A	312	13
	SoundRecover <i>non-linear frequency compression</i>	■	■	■
	Programmable Channels	6	6	6
	SoundFlow Standard <i>adapts for seamless transition</i>	■	■	■
Standard Instrument Specific Features <ul style="list-style-type: none"> DataLogging <i>(records and stores personal use and listening preferences)</i> End of battery life indicator 	digital AudioZoom <i>automatic directional microphone</i>		■	■
	WhistleBlock <i>more stable amplification without sound quality impact</i>	■	■	■
	BassBoost <i>increased low frequency gain and output</i>	■ <i>(CIC/MC P only)</i>	■ <i>(ITC/HS P only)</i>	■ <i>(FS P only)</i>
	NoiseBlock <i>detects and reduces annoying sounds</i>	■ <i>(6 Channels)</i>	■ <i>(6 Channels)</i>	■ <i>(6 Channels)</i>
	AOV - Acoustically Optimized Vent <i>minimizes occlusion with maximum vent size</i>	■	■	■
	EasyPhone - Automatic Telephone		■ <i>(Optional)</i>	■ <i>(Optional)</i>
	Telecoil		■	■
	Number of Manual Programs	2	2+2 dedicated T	2+2 dedicated T
Options <ul style="list-style-type: none"> Choice of BTE colors:  Choice of micro BTE colors:  Choice of micro Petite BTE colors:  Choice of Custom colors: See page 8 	EasyAudio - <i>binaural stereo streaming</i> EasyBluetooth™ - <i>handsfree cell phone use</i>		■ <i>(Wireless only)</i>	■
	QuickSync <i>binaural volume and program control</i>		■ <i>(Wireless only)</i>	■
	iCom - <i>wireless Bluetooth™ Et FM</i>		■ <i>(Wireless only)</i>	■
	iCube <i>CableFree Fitting</i>		■ <i>(Wireless only)</i>	■
	myPilot (remote control) <i>Bi-directional access</i>		■ <i>(Wireless only)</i>	■
	KeyPilot2 Et WatchPilot2 <i>(remote controls)</i>	■	■	■
	Design Integrated FM Receiver			
	Universal FM Receiver <i>myLink</i>		■	■
	FM Transmitters <i>inspiro, SmartLink+, EasyLink+, ZoomLink+</i>		■	■
	Tamper proof Battery compartment			


Certéna Art



Features	BTE M	BTE P	BTE SP	Micro
	Wireless	Wireless	Wireless	Wireless/Petite
Fitting Range	53 dB/125 dB SPL	63 dB/130 dB SPL	70 dB/131 dB SPL	46 dB/123 dB SPL - Petite 53 dB/124 dB SPL - Wireless
Battery Size	13	13	13	10A - Petite 13 - Wireless
SoundRecover <i>non-linear frequency compression</i>	■	■	■	■
Programmable Channels	6	6	6	6
SoundFlow Standard <i>adapts for seamless transition</i>	■	■	■	■
digital AudioZoom <i>automatic directional microphone</i>	■	■	■	■
WhistleBlock <i>more stable amplification without sound quality impact</i>	■	■	■	■
BassBoost <i>increased low frequency gain and output</i>			■	
NoiseBlock <i>detects and reduces annoying sounds</i>	■ (6 Channels)	■ (6 Channels)	■ (6 Channels)	■ (6 Channels)
AOV - Acoustically Optimized Vent <i>minimizes occlusion with maximum vent size</i>				
EasyPhone - Automatic Telephone	■	■	■	■ (Wireless only)
Telecoil	■	■	■	■
Number of Manual Programs	2+2 dedicated T/FM	2+2 dedicated T/FM	2+2 dedicated T/FM	2+2 dedicated T - Petite 2+2 dedicated T/FM - Wireless
EasyAudio - binaural stereo streaming EasyBluetooth™ - handsfree cell phone use	■	■	■	■ (Wireless only)
QuickSync <i>binaural volume and program control</i>	■	■	■	■ (Wireless only)
iCom - wireless Bluetooth™ Et FM	■	■	■	■ (Wireless only)
iCube <i>CableFree Fitting</i>	■	■	■	■ (Wireless only)
myPilot (remote control) <i>Bi-directional access</i>	■	■	■	■ (Wireless only)
KeyPilot2 Et WatchPilot2 (remote controls)	■	■	■	■
Design Integrated FM Receiver	ML9i	ML9i	ML9i	ML12i (Wireless only)
Universal FM Receiver <i>myLink</i>	■	■	■	■
FM Transmitters <i>inspiro, SmartLink+, EasyLink+, ZoomLink+</i>	■	■	■	■
Tamper proof Battery compartment	■ (Optional)	■ (Optional)	■ (Optional)	■ (Wireless option only)

Milo Plus



	Features	CIC/MC	ITC/HS dAZ	FS dAZ FS P dAZ
Fitting Recommendations <ul style="list-style-type: none"> Price-sensitive clients Easy to fit Accommodates virtually all degrees of hearing loss 	Fitting Range	35dB/113dB SPL	48 dB/113 dB SPL	50 dB/115 dB SPL - FS 55 dB/122 dB SPL - FS P
	Battery Size	10A	312	13
	Programmable Channels <i>4 Channel Gain and MPO Shaping</i>	4	4	4
	Directionality <i>dual microphone system</i>		■	■
Standard Instrument Specific Features <ul style="list-style-type: none"> Fully programmable digital hearing system DataLogging Mute mode & Delayed Start-up End of battery-life indicator Audiogram Direct 	WhistleBlock <i>more stable amplification without sound quality impact</i>	■	■	■
	BassBoost <i>increased low frequency gain and output</i>			■ <i>(FS P Only)</i>
	Wind & weather protector			
	NoiseBlock <i>detects and reduces annoying sounds</i>	■ <i>(4 Channels)</i>	■ <i>(4 Channels)</i>	■ <i>(4 Channels)</i>
	AOV - Acoustically Optimized Vent <i>minimizes occlusion with maximum vent size</i>	■	■	■
	Telecoil		■	■
	Number of Manual Programs	2	2+2 dedicated T	2+2 dedicated T
Options <ul style="list-style-type: none"> Choice of BTE colors:  Choice of Custom colors: See page 8 	Volume control	■ <i>(Optional)</i>	■ <i>(Optional)</i>	■ <i>(Optional)</i>
	Acoustic confirmation for all manual changes	■	■	■
	Program Switch	■ <i>(Optional)</i>	■ <i>(Optional)</i>	■ <i>(Optional)</i>
	Design Integrated FM Receiver			
	Universal FM Receiver <i>myLink+</i>		■	■
	FM Transmitters <i>inspiro, SmartLink+, EasyLink+, ZoomLink+</i>		■	■
	Tamper proof battery compartment			




Features	Milo Plus micro	Milo Plus SP	Milo Plus UP
Fitting Range	50 dB/127 dB SPL	65 dB/133 dB SPL	75 dB/138 dB SPL
Battery Size	13	13	675
Programmable Channels <i>4 Channel Gain and MPO Shaping</i>	4	4	4
Directionality <i>dual microphone system</i>	■	■	■
WhistleBlock <i>more stable amplification without sound quality impact</i>	■	■	■
BassBoost <i>increased low frequency gain and output</i>		■	■
Wind & weather protector	■	■	■
NoiseBlock <i>detects and reduces annoying sounds</i>	■ <i>(4 Channels)</i>	■ <i>(4 Channels)</i>	■ <i>(4 Channels)</i>
AOV - Acoustically Optimized Vent <i>minimizes occlusion with maximum vent size</i>			
Telecoil	■	■	■
Number of Manual Programs	2+2 dedicated T/FM	2+2 dedicated T/FM	2+2 dedicated T/FM
Volume control		■	■
Acoustic confirmation for all manual changes	■	■	■
Program Switch	■	■	■
Design Integrated FM Receiver	ML12i	ML11i	ML10i
Universal FM Receiver <i>MLxi and myLink+</i>	■	■	■
FM Transmitters <i>inspiro, SmartLink+, EasyLink+, ZoomLink+</i>	■	■	■
Tamper proof battery compartment	■ <i>(Optional)</i>	■ <i>(Optional)</i>	■ <i>(Optional)</i>

Audéo ZIP



	Features	IX	V	III
	Fitting Range	40 dB/112 dB SPL	40 dB/112 dB SPL	40 dB/112 dB SPL
	Battery Size	10A	10A	10A
	SoundRecover <i>non-linear frequency compression</i>	■	■	■
Fitting Recommendations	Programmable Channels	20	16	6
<ul style="list-style-type: none"> ■ Ideal for first time users ■ A nearly invisible hearing solution for mild-to-moderate hearing loss 	SoundFlow <i>adapts for seamless transition</i>	Premium	Advanced	Standard
	Directionality	Omni	Omni	Omni
	WhistleBlock <i>more stable amplification without sound quality impact</i>	■	■	■
	EchoBlock <i>reverberation canceler</i>	■		
Standard Instrument Specific Features	WindBlock <i>suppresses wind noise while maintaining directionality</i>	■	■	
<ul style="list-style-type: none"> ■ InstantFit CIC ■ SelfLearning and DataLogging with SmartTips (<i>Records and stores personal use and listening preferences</i>) (<i>DataLogging only - Audéo ZIP III</i>) ■ End of battery life indicator 	NoiseBlock <i>detects and reduces annoying sounds</i>	■ (20 Channels)	■ (16 Channels)	■ (6 Channels)
	Number of Manual Programs <i>accessible via remote only</i>	3	3	2
	SoundRelax <i>decreases unwanted sharp sounds</i>	■		
	TwinVent <i>minimizes occlusion 5mm vent size equivalent</i>	■	■	■
	Closed/Open/Power Dome Fitting	■	■	■
<ul style="list-style-type: none"> ■ 1 Dome ■ 2 Receiver housing with left/right indicator ■ 3 MultiFlex Canal ■ 4 SecureFit Tab ■ 5 Microphone with protector ■ 6 Removal pin ■ 7 Battery compartment ■ 8 Serial number (<i>inside battery compartment</i>) 	KeyPilot2 Et WatchPilot2 (remote controls)	■	■	■

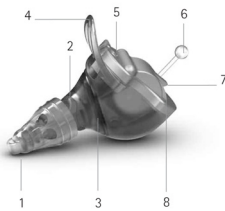

- Choice of Faceplate colors:
See page 8

Custom Colors

Faceplate



Clear 21

Shell



Audéo ZIP InstantFit CIC Colors

