

Phonak Audéo™ P

Product information



Phonak Audéo P is the comprehensive RIC portfolio for mild to profound hearing losses. Phonak Audéo P includes four designs, four external receivers and four performance levels. Powered by the innovative technology of the Phonak Paradise platform, Phonak Audéo P hearing aids are the ideal solution for those clients that prefer best hearing performance in flexible and discreet RIC housings.

For more information visit www.phonakpro.com.au

Phonak Audéo P charging options



Phonak Charger Case Combi

The Phonak Charger Case Combi is a charger, drying kit and protective hard case all in one.



Phonak Power Pack

To charge rechargeable hearing aids on the go, the Phonak Power Pack can be easily attached to the bottom of the Phonak Charger Case Combi.



Phonak Mini Charger Case RIC

The Phonak Mini Charger Case, with standard USB-C charging outlet, is a compact charger for two hearing aids.



Phonak Charger BTE RIC

This charger is suitable for all Paradise rechargeable BTE and RIC hearing aids. It features a holding space for large earmoulds, a USB-C charging inlet and an indicator light.

Phonak Audéo P Bluetooth® versions

Direct calls/audio streaming

Bluetooth® 4.2 wireless technology and most older Bluetooth phones

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.

Wireless communication portfolio

	Audéo P-312 Audéo P-312 Trial P90/P70/P50/P30	Audéo P-13T Audéo P-13T Trial P90/P70/P50/P30	Audéo P-R Audéo P-R Trial P90/P70/P50/P30	Audéo P-RT Audéo P-RT Trial P90/P70/P50/P30
Wireless accessory				
TV Connector	•	•	•	•
Phonak PartnerMic	•	•	•	•
Phonak RemoteControl	•	•	•	•
Roger				
RogerDirect™	•	•	•	•
Apps				
myPhonak	•	•	•	•

For more information please visit www.phonakpro.com/audeoparadise

For information on Roger microphones please visit www.phonakpro.com/roger

For more information please visit www.phonakpro.com/smartapps

Product description

					
		Audéo P-312 Audéo P-312 Trial P90/P70/P50/P30	Audéo P-13T Audéo P-13T Trial P90/P70/P50/P30	Audéo P-R Audéo P-R Trial P90/P70/P50/P30	Audéo P-RT Audéo P-RT Trial P90/P70/P50/P30
Battery size		312	13		
Multi-function button		•	•	•	•
Telecoil			•		•
Bluetooth®		•	•	•	•
Nano coating		•	•	•	•
IP rating		IP68 ¹	IP68 ¹	IP68 ¹	IP68 ¹
Dimensions	L x W x D	26.8 x 11.2 x 7.1 mm (1.05 x 0.44 x 0.28")	29.0 x 11.9 x 8.1 mm (1.13 x 0.47 x 0.32")	28.0 x 11.9 x 8.0 mm (1.11 x 0.47 x 0.32")	30.5 x 12.5 x 8.0 mm (1.20 x 0.49 x 0.32")
	Weight	1.1 g (0.04 oz)	1.4 g (0.05 oz)	2.1 g (0.07 oz)	2.4 g (0.08 oz)
Charging times ²	0% to 100%			3 h	3 h
	0% to 80%			1.5 h	1.5 h
	30 min. fitting			10 min.	10 min.
S Receiver					
Max. Power Output (dB SPL)	2 cc coupler	111	111	111	111
Max. gain (dB)	2 cc coupler	47	47	47	47
Frequency range (Hz)		<100 - >8000	<100 - >8000	<100 - >8000	<100 - >8000
Working current (mA)		1.4	1.4		
M Receiver					
Maximum Power Output (dB SPL)	2 cc coupler	114	114	114	114
Max. gain (dB)	2 cc coupler	51	51	51	51
Frequency range (Hz)		<100 - >8000	<100 - >8000	<100 - >8000	<100 - >8000
Working current (mA)		1.4	1.4		
P Receiver					
Maximum Power Output (dB SPL)	2 cc coupler	122	122	122	122
Max. gain (dB)	2 cc coupler	59	59	59	59
Frequency range (Hz)		<100 - 6300	<100 - 6300	<100 - 6300	<100 - 6300
Working current (mA)		1.5	1.5		
UP Receiver					
Max. Power Output (dB SPL)	2 cc coupler	132	132	132	132
Max. gain (dB)	2 cc coupler	71	71	71	71
Frequency range (Hz)		<100 - 6100	<100 - 6100	<100 - 6100	<100 - 6100
Working current (mA)		1.4	1.4		

Housing colours

						
P1 Sand Beige	H0 Beige	P4 Chestnut	P5 Champagne	P6 Silver Gray	P7 Graphite Gray	P8 Velvet Black

1 IP68 indicates that the hearing aid is water resistant and dust tight. It survived continuous immersion in 1 metre of water for 60 minutes and 8 hours in a dust chamber as per the IEC60529 standard.

2 Actual charging time depends on the remaining battery life, but it will not exceed 3 hours.

Performance levels

	Premium P90	Advanced P70	Standard P50	Essential P30
AutoSense OS™ 4.0	Premium	Advanced	Standard	Essential
Comfort in Echo	•			
Speech in Car	•			
Speech in Loud Noise	•			
Music	•	•		
Comfort in Noise	•	•	•	
Speech in Noise	•	•	•	•
Calm Situation	•	•	•	•
Media speech + mic	•	•	•	•
Media music + mic	•	•	•	•
Additional programs				
Max. additional programs	4	4	4	4
Comfort in Echo	•			
Speech in Loud Noise	•	•		
Speech in 360°	•	•		
Speech in Noise	•	•	•	•
Calm Situation	•	•	•	•
Comfort in Noise	•	•	•	•
Music	•	•	•	•
Acoustic phone	•	•	•	•
Custom program	•	•	•	•
Streaming programs				
RogerDirect + mic	•	•	•	•
PartnerMic + mic	•	•	•	•
Phone call + mic	•	•	•	•

	Premium P90	Advanced P70	Standard P50	Essential P30
Features	Premium	Advanced	Standard	Essential
UltraZoom				
Dynamic Noise Cancellation	•			
Speech Enhancer	•			
EchoBlock	•			
Tap control*				
- Siri/Google Assistant	•	•		
- Pause/resume streaming	•	•		
- Accept/end phone call	•	•		
StereoZoom	•	•		
WindBlock	•	•		
SoundRelax	•	•		
Motion Sensor Hearing*	•	•	•	
DuoPhone	•	•	•	
Real Ear Sound	•	•	•	
Roger and Directional setting	•	•	•	
SoundRecover2	•	•	•	•
WhistleBlock	•	•	•	•
NoiseBlock	•	•	•	•
RogerDirect™	•	•	•	•
User Preference Tuning	•	•	•	•
QuickSync	•	•	•	•
AOV	•	•	•	•
Tinnitus Balance	•	•	•	•
auto Acclimatisation	•	•	•	•
Environmental balance	•	•	•	•
Fine tuning channels	20	20	16	12

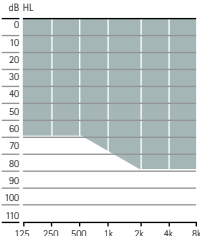
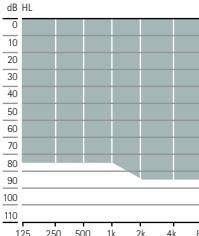
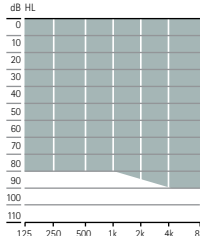
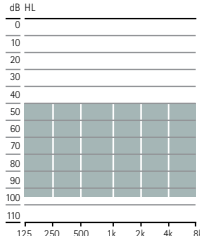
* Only available for Phonak Audéo P-R and RT

Fitting

	Audéo P-312 Audéo P-312 Trial P90/P70/P50/P30	Audéo P-13T Audéo P-13T Trial P90/P70/P50/P30	Audéo P-R Audéo P-R Trial P90/P70/P50/P30	Audéo P-RT Audéo P-RT Trial P90/P70/P50/P30
Phonak Target	7.0 or higher	7.0 or higher	7.0 or higher	6.1 or higher
Noahlink Wireless	•	•	•	•
NOAHLink		•		
HI-PRO, HI-PRO2		•		

The Bluetooth® word mark and logos are registered trademarks owned by the Bluetooth SIG, Inc.
IOS is a trademark of Cisco, Inc.
Android is a trademark of Google, Inc.

Acoustic coupling

Length	S Receiver	M Receiver	P Receiver	UP Receiver
Available in four different lengths (0, 1, 2, 3), left and right	•	•	•	•
Dome				
Cap dome	•	•	•	
Open dome Available in three sizes (S, M, L)	•	•	•	
Vented dome Available in three sizes (S, M, L)	•	•	•	
Power dome Available in three sizes (S, M, L)	•	•	•	
SlimTip*				
Acryl	•	•	•	
Silicone	•	•	•	
cShell				
Acryl	•	•	•	•
Titanium	•	•	•	•
Fitting ranges				
				
	Mild to moderately-severe hearing loss.	Mild to severe hearing loss.	Mild to severe hearing loss.	Moderate to profound hearing loss.

Phonak Audéo P Trial™

The new Phonak Audéo P Trial hearing aids can be quickly set to the required performance level in order to suit your recommendation. All levels of technology in one Phonak Audéo P Trial hearing aid provides you with many options for providing the hearing solution at the first appointment or bridging the gap during repairs.

*Available with extra canal length option

life is on

At Phonak, we believe that well-hearing equates to well-being and is essential to living life to the fullest. For more than 70 years, we have remained passionate about creating a world where 'life is on' for everyone. Our innovative hearing solutions are designed for people of all ages and all degrees of hearing loss, to connect socially, thrive mentally and emotionally.

www.phonakpro.com.au



027-0602-02/V1.00/2020-04/na Printed in Australia © 2020 Sonova AG All rights reserved