


Phonak CROS II

The smart solution for single-sided deafness



Introducing Phonak CROS II

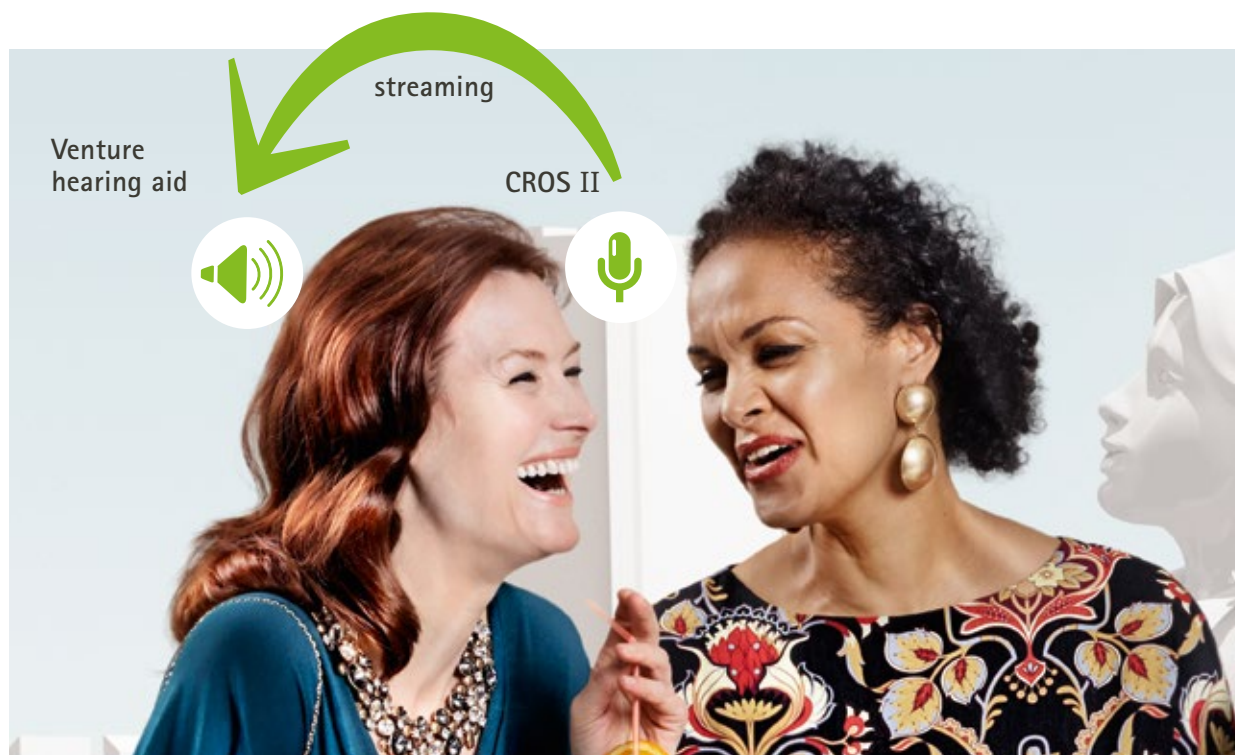

Phonak CROS II is the discreet wireless CROS transmitter based on the Venture platform. Featuring the unique Phonak Binaural VoiceStream Technology™, it wirelessly transmits the sound from the unaidable ear to a Venture hearing aid on the better hearing ear. As a result a person with single-sided deafness (SSD) can hear speech from the side they have not heard from before.

Phonak CROS II is easier to fit, delivers better speech understanding and consumes up to 30% less battery than the previous Phonak CROS.



Phonak
CROS II-312

Phonak
CROS II-13



Phonak
CROS II-312
Custom



Phonak
CROS II-13
Custom

Easy to fit

We recognize that you may not be performing a CROS fitting weekly. But do not worry, the streamlined fitting process of a Phonak CROS II system is as easy and intuitive as a binaural fit.



Go to phonakpro.com/cros to see how simple it is to fit.



iCube II fitting device

Easy to demonstrate


When using iCube II, Phonak CROS II is active during the fitting. Your clients can **instantly experience the benefits** of the CROS solution and provide feedback during the fitting session.



CROS II Hook for CROS II-13 only

Easy to secure

The CROS II can be fit in many configurations to meet the needs of your client. From the security of a tonehook-earmold combination to the sleek yet secure design of the custom CROS Tip, Phonak has the solution.

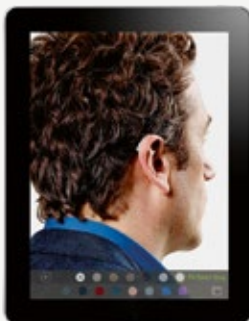

Show your clients what can't be seen

With the Virtual Mirror App you can use an iPad to mirror your client's ear with Phonak CROS II. Virtually change the color according to their preference and show them high definition images.

Download the Phonak Virtual Mirror App for free from the App Store.



More speech understanding

For the first time ever, **people with SSD can easily follow conversations even in** the most challenging situation of all: **noisy surroundings.**

This is thanks to StereoZoom, the binaural beamforming feature that can focus on a single voice in a crowd and reduce background noise. As a result, the CROS II transmitter streams less noise to the hearing aid on the other side, improving speech understanding in noisy situations.



Go to phonakpro.com/cros to hear the superior performance of CROS II in noisy surroundings.

Longer battery life

CROS II constantly transmits audio signals from the unaidable ear to the hearing aid on the better ear. Thanks to the Venture chip with more efficient wireless processing, your clients will profit from a **30% longer battery life than the previous Phonak CROS.**

CROS volume control

Some clients prefer to be in control and be able to **separately adjust the volume with the CROS transmitter.** This is possible both with the on-board push button or via an App on your Smartphone.




With the RemoteControl App you can control the hearing aid and the CROS volume using your Smartphone.

Download the Phonak RemoteControl App for free from the App Store or from phonakpro.com/rcapp.

One CROS II, many options

The performance of Phonak CROS II in noisy environments improves with the performance level of the hearing aid it works with. Combine with a V70 or a V90 to **benefit from the best performance possible in noise, thanks to StereoZoom.**

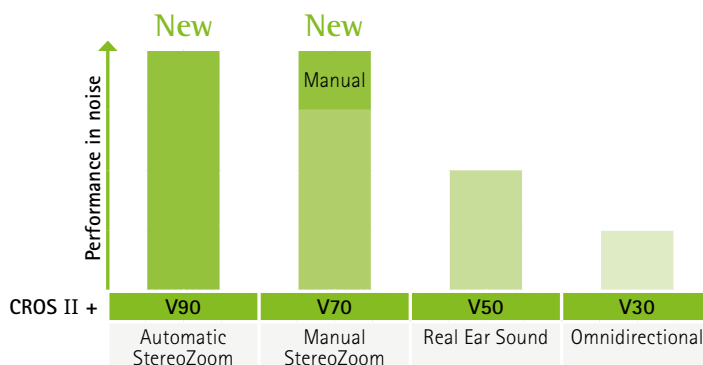

Combine CROS II with any Phonak Venture hearing aid and choose from the same color range as the Venture hearing aid family.

For a list of all shell and faceplate colors for the CROS II Custom products, see the CROS II Product Information.



Life is on

We are sensitive to the needs of everyone who depends on our knowledge, ideas and care. And by creatively challenging the limits of technology, we develop innovations that help people hear, understand and experience more of life's rich soundscapes.

**Interact freely. Communicate with confidence.
Live without limit. Life is on.**

www.phonakpro.com

