


Phonak TVLink basestation

Technical Data

Phonak TVLink basestation, the Bluetooth® transmitter that turns Phonak hearing instruments into the lightest wireless TV headset ever. The TVLink basestation wirelessly sends audio signals from the TV or any other audio source to iCom within a range of up to 30 meters (100 ft). The unique StereoSound guarantees a perfect listening experience and the amazingly short transmission delay ensures listening pleasure even for open fittings. The TVLink basestation also serves as a charging cradle for iCom.



Description

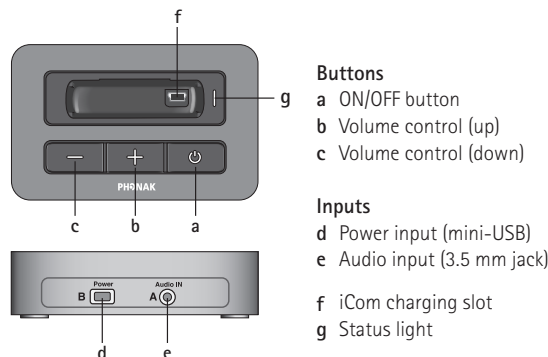

Special features

- Compatible with all wireless Phonak hearing instruments
- Easy and quick installation
- Complete set with power plug and cables
- No line of sight to hearing instruments necessary
- Up to 30 meters (100 ft) transmission range
- Digital volume adjustment with audio-visual feedback
- Embedded charging slot for iCom
- Power saving option

General information

Dimensions:	L x W x H: 90 x 60 x 23 mm (3.8 x 2.5 x 1 inches)
Weight:	100 g
Color:	Silver metallic
Operating conditions:	The product has been designed for trouble-free operation without restrictions when used as intended, unless otherwise stated in the user guide -10°/55° Celsius and relative humidity of <95% (non condensing)
Transport and storage conditions:	-20°/60° Celsius and relative humidity of 90% for a long period of time

Status light indication

- Power ON
- Transmitting to iCom
- - - Pause pressed on iCom
- Bluetooth pairing

Manual Bluetooth pairing:
- Press ON/OFF button for 5 seconds

Delete Bluetooth pairing table:
- Press ON/OFF button for 10 seconds



Technical Data

Phonak TVLink basestation

Power supply

Type:	DC, stabilized
Secondary voltage:	5 VDC
Max. current:	500 mA
Plug:	Mini-USB
Important notice:	Only use the original Phonak power plug or an USB certified power adaptor > 500 mA

Audio information

Audio input:	Wired (3.5 mm socket)
Audio quality:	Stereo, 20 Hz – 22 kHz

Based on Bluetooth® technology

Bluetooth profile:	
A2DP, Advanced Audio Distribution Profile	<ul style="list-style-type: none">▪ Optimized for use with Phonak iCom▪ Compliance to non-Phonak devices is not guaranteed
Frequency range:	2402.0 to 2480.0 MHz
RF power:	100 mW max.
Bluetooth range:	Up to 30 m (100 ft)
Bluetooth compliance:	Bluetooth version 2.1, Class 1
Maximum paired devices:	4
QDID:	B016567

Standards

Europe

Radio:	EN 300 328
EMC:	EN 301.489
El. Safety:	EN 60950-1

USA	CFR 47, part 15
-----	-----------------

Canada	RSS-210 Issue 7
--------	-----------------



The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Phonak is under license. Other trademarks and trade names are those of their respective owners.



This device is a class 2 device and can be operated in the EU member states without restrictions for indoor use, but must not operated outdoors in France.

Accessories

- Audio cable 3.5 mm plug
- Audio cable RCA/Cinch
- SCART adaptor (EU only)
- Power plug: country specific
- User guide
- Quick setup guide